


# AMRAD MEDICAL

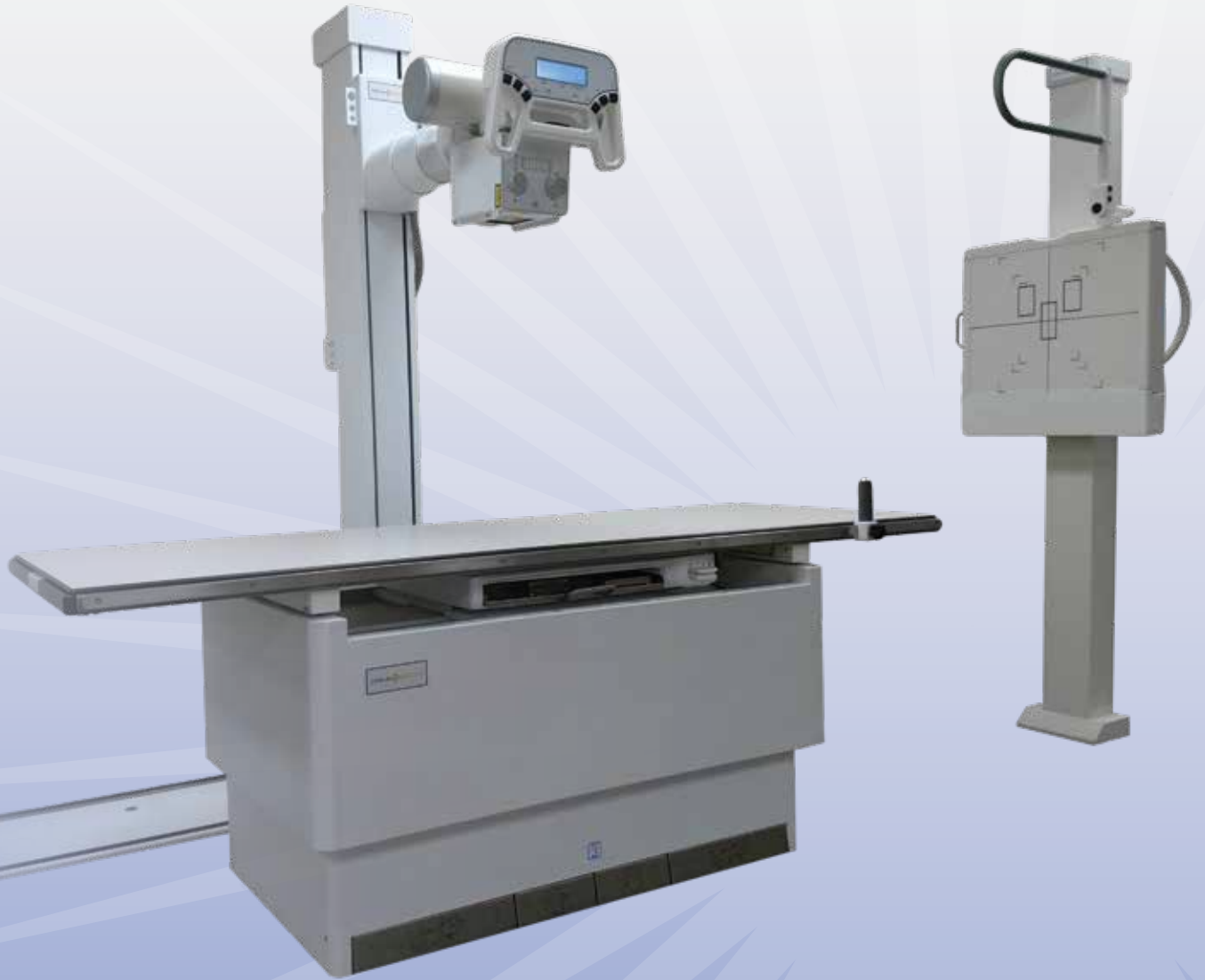

**DFMT | FMT**

Premier Radiographic Systems

[www.AmRadMedical.com](http://www.AmRadMedical.com)

# Features & Benefits




DFMT



Fast and Accurate Patient Positioning

## Digital Floor Mounted Tubestand

- Easy to read LCD shows SID continuous to table-top and table receptor with discrete SID of 40" and 72" to the wallstand
- Cross-table views are easily accomplished using the rotational and transverse positioning capabilities
- 10 ft. floor track allows for off-table and weight bearing exams and is customizable down to 8 ft.

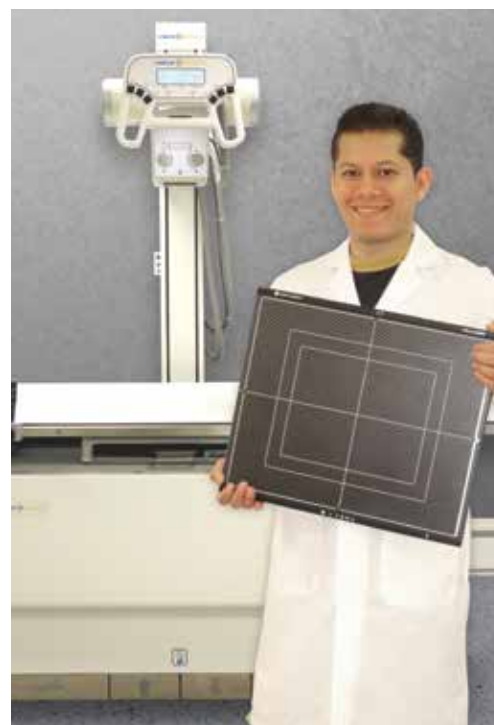

Compatible with most DR Solutions



Easy Patient Transfer



Convenient Lock-Release Pedals



Standing Knees



Ergonomic Design

## Elevating Four-way Float Top Table

- Flat, four-way, floating top permits easy transfer and positioning of patients
- Designed for all types of applications, available in 400lb and 650lb patient capacities
- Four-pedal design allows for easy accessibility to table lock release
- Collision-protection and fail-safe locks ensure safe operation
- Four-way lock release and elevation handle (DFMT)

### Options include:

- 35.5" wide table top
- Four-way lock release handle
- Non-elevating version available in 400 lb and 650 lb patient capacities
- Rotating Cassette Tray

## Wallstand

- Center mounting or side mounting of Bucky makes this wallstand flexible for any patient application or installation limitation
- Ergonomically designed hand control for ease of positioning wall receptor
- Side-mounted handgrips stabilize patients during PA exams
- Overhead Grip (DFMT)
- Remote Collimator Light

### Options include:

- Dual-sided Lock Release
- Rotating Cassette Tray



35.5" Wide Table Top



Economical Alternative FMT



Table Lock Release Handle with Elevation (FMT, optional)



Overhead Grip for J1000 Wallstand (FMT, optional)

## System Options (L to R, above)

The AmRad™ Medical DFMT and FMT Systems can be configured to meet the needs of any department. Options include, but are not limited to:

- An extra-wide table top for the S222 elevating table to help patients feel more secure
- Four-way lock release and elevation handle for easy patient positioning can be located anywhere on the rails of the S223 non-elevating or S222 elevating tables (DFMT, standard; FMT, optional)
- The FMT is an economical alternative to the DFMT with the same versatility and functionality as the DFMT
- Patient Overhead Grip for use with the J1000 wallstand helps steady patients for upright exams (DFMT, standard; FMT, optional)



SAI GLOBAL  
ISO 13485  
Quality

**AMRAD**  **MEDICAL**

[www.AmRadMedical.com](http://www.AmRadMedical.com)